**medi GmbH & Co. KG**  
**Medicusstraße 1**  
**95448 Bayreuth / Germany**  
**T +49 921 912-0**  
**F +49 921 912-780**  
**ortho@medi.de**  
**www.medi.de/en**

medi Australia Pty Ltd  
 83 Fennell Street  
 North Parramatta NSW 2151  
 Australia  
 T +61-2 9890 8696  
 F +61-2 9890 8439  
 sales@medi australia.com.au  
 www.medi australia.com.au

medi Austria GmbH  
 Adamgasse 16/7  
 6020 Innsbruck  
 Austria  
 T +43 512 57 95 15  
 F +43 512 57 95 15 45  
 vertrieb@medi-austria.at  
 www.medi-austria.at

medi Belgium NV  
 Staatsbaan 77/0099  
 3945 Ham  
 Belgium  
 T: + 32-11 24 25 60  
 F: +32-11 24 25 64  
 info@medibelgium.be  
 www.medibelgium.be

medi Canada Inc / médi Canada Inc  
 104-1375 Lionel-Boulet,  
 Varennes, Québec,  
 QC Canada J3X 1P7  
 T +1 450-583-3317 / +1 800 361 3153  
 F +1 888-583-6827  
 service@medicanada.ca  
 www.medicanada.ca

MAXIS a.s.,  
 medi group company  
 Slezská 2127/13  
 120 00 Prague 2  
 Czech Republic  
 T: +420 571 633 510  
 F: +420 571 616 271  
 info@maxis-medica.com  
 www.maxis-medica.com

medi Danmark ApS  
 Vejlegardsvej 59  
 2665 Vallensbaek Strand  
 Denmark  
 T +45-70 25 56 10  
 kundeservice@sw.dk  
 www.medit danmark.dk

medi Bayreuth Espana SL  
 C/Canigo 2-6 bajos  
 Hospitalet de Llobregat  
 08901 Barcelona  
 Spain  
 T +34-932 60 04 00  
 F +34-932 60 23 14  
 medi@mediespana.com  
 www.mediespana.com

medi France  
 Z.I. Charles de Gaulle  
 25, rue Henri Farman  
 93297 Tremblay en France Cedex  
 France  
 T +33-1 48 61 76 10  
 F +33-1 49 63 33 05  
 infos@medi-france.com  
 www.medi-france.com

medi Hungary Kft.  
 Bokor u. 21.  
 1037 Budapest  
 Hungary  
 T +36 1 371-0090  
 F +36 1 371-0091  
 info@medi.hu  
 www.medi.hu

medi Japan K.K.  
 5-1 Nihonbashi Kabutocho Chuoku  
 103-0026 Tokyo  
 Japan  
 T: +81 3 6778 2590  
 F: +81 3 5847 7901  
 info@medi-japan.jp  
 www.medi-japan.co.jp

medi Nederland BV  
 Heusing 5  
 4817 ZB Breda  
 The Netherlands  
 T +31-76 57 22 555  
 F +31-76 57 22 565  
 info@medi.nl  
 www.medi.nl

medi Polska Sp. z o.o.  
 ul. Łabędzka 22  
 44-121 Gliwice  
 Poland  
 T: +48-32 230 60 21  
 F: +48-32 202 87 56  
 info@medi-polska.pl  
 www.medi-polska.pl

medi Bayreuth Unipessoal, Lda  
 Rua do Centro Cultural, no. 43  
 1700-106 Lisbon  
 Portugal  
 T +351-21 843 71 60  
 F +351-21 847 08 33  
 medi.portugal@medibayreuth.pt  
 www.medi.pt

medi RUS LLC  
 Business Center NEO GEO  
 Butlerova Street 17  
 117342 Moscow  
 Russia  
 T +7-495 374 04 56  
 F +7-495 374 04 56  
 info@medius.ru  
 www.medi us.ru

medi Sweden AB  
 Box 6034  
 192 06 Sollentuna  
 Sweden  
 T +46 8 96 97 98  
 F +46 8 626 68 70  
 info@medi.se  
 www.medi.se

medi Trading (Shanghai) Co. Ltd.  
 Room 1703, T1, CES WEST BUND CENTER,  
 No. 277, Long Lan Road  
 Xuhui District  
 200232 Shanghai  
 People's Republic of China  
 T: +86-21 50582319  
 F: +86-21 50582319


Medi Turk Ortopedi Medikal Iç ve Dis Tic. A.Ş.  
 Mustafa Kemal Mah. 2156 Sokak No: 16/4-5  
 06510 Çankaya Ankara  
 Turkey  
 T: +90 312 435 20 26  
 F: +90 312 434 17 67  
 info@medi-turk.com  
 www.medi-turk.com

medi Ukraine LLC  
 Evhena Sverstiuka str, 11.  
 Kiev 02002  
 Ukraine  
 T: +380 44 591 11 63  
 F: +380 44 392 73 73  
 info@medi.ua  
 http://medi.ua

medi USA L.P.  
 6481 Franz Warner Parkway  
 Whitsett, N.C. 27377-3000  
 USA  
 T +1-336 4 49 44 40  
 F +1-888 5 70 45 54  
 info@mediusa.com  
 www.mediusa.com

UK Responsible Person  
 medi UK Ltd.  
 Plough Lane  
 Hereford HR4 0EL  
 Great Britain  
 T +44-1432 37 35 00  
 F +44-1432 37 35 10  
 enquiries@mediuk.co.uk  
 www.mediuk.co.uk

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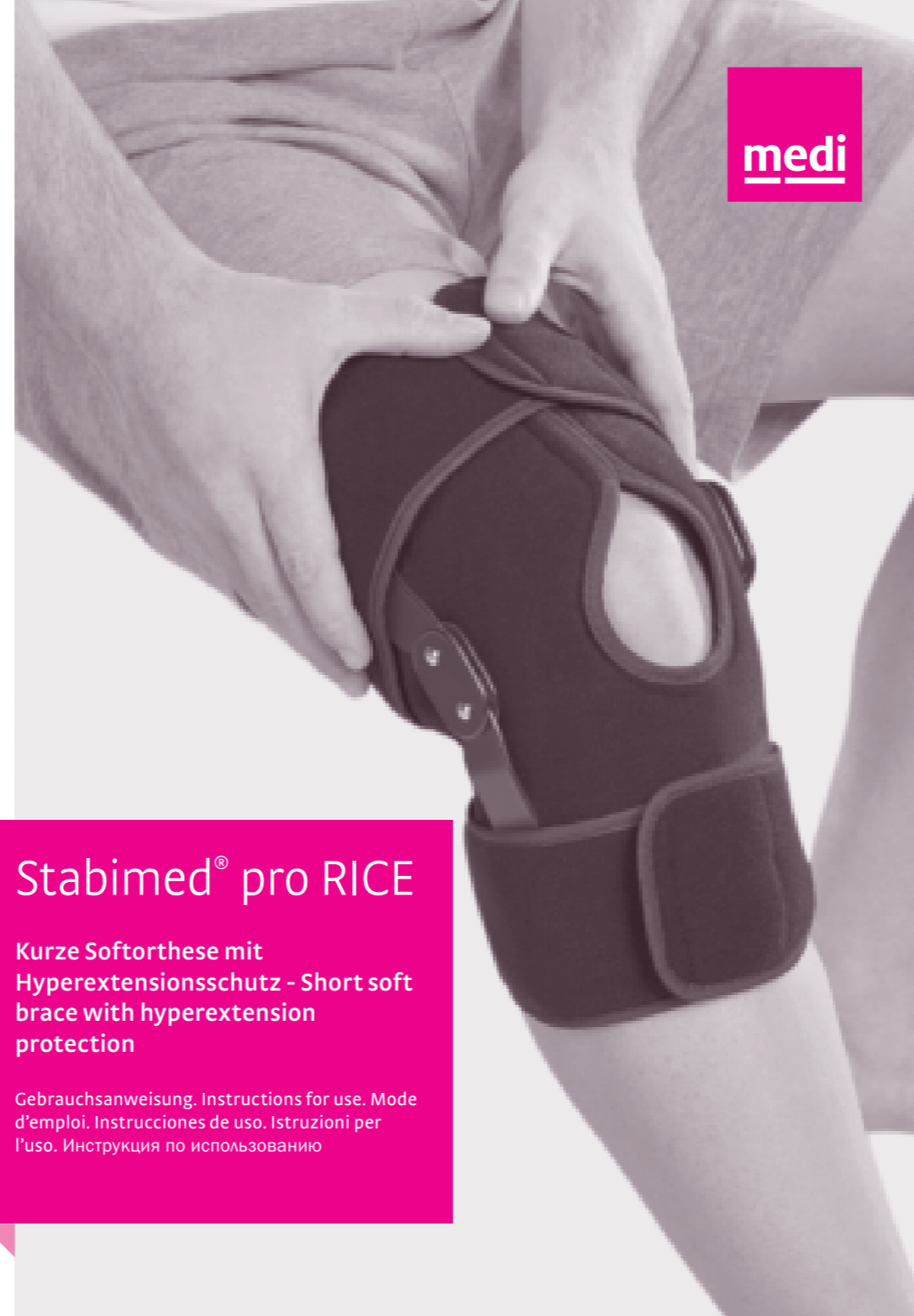### Important notes

This medical device **MD** is made only for multiple uses on a single patient. If it is used for treating more than one patient, the manufacturer's product liability will become invalid. If you experience excessive pain or any discomfort whilst wearing, remove the product and immediately contact your doctor or medical supply store. Only wear the product on unbroken skin or over injured skin that has been medically treated, not directly on injured or damaged skin, and only after previous consultation with a medical professional.

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# Stabimed® pro RICE

## Intended purpose

The Stabimed pro RICE is a knee support brace without limitation of extension / flexion.

## Indications

All indications that require reinforcement of the physiological stabilisation of the knee joint such as:

- Knee instabilities
- Following meniscus injuries
- Following endoprosthetic intervention
- Osteoarthritis of the knee

## Contraindications

None known at present.

## Risks / Side Effects

If the aids are too tight, this can cause local pressure sores or constriction of blood vessels or nerves. Therefore you should consult the doctor treating you before you use it if the following circumstances apply to you:

- Condition or injury of the skin in the area of application, particularly with infectious indications (excessive heat, swelling or reddening)
- Paraesthesia and circulatory disorders (e.g. as seen in diabetes, varicose veins)
- Lymph drainage disorders – equally, swelling of soft tissue adjacent to the area of application which is not related to a condition

When wearing snug-fitting aids, this can result in localised skin inflammation or irritation that is due to mechanical irritation of the skin (particularly in conjunction with sweating) or due to the composition of the material.

## Intended users and patient target groups

Healthcare professionals and patients, including persons who play a supporting role in care, are the intended users for this device, provided that they have received briefing on this from healthcare professionals.

Patient target group: Healthcare professionals should provide care to adults and children, applying the available information on the measurements/sizes and required functions/indications, in line with the information provided by the manufacturer, and acting under their own responsibility.

## Important notes

Please ensure that you only use this brace without its articulated joints if the doctor treating you has told you to do so.

## Fitting instructions

- Open both straps (Fig. 1). The valve is at the top right and the straps are also facing right.
- Bend your leg to approximately 45°. Place the brace onto the front of the leg so that the lower edge of the hole for the kneecap is sitting on the lower side of the kneecap and is positioned centrally (Fig. 2).
- Now position both of the joints on the brace. To do this, put the joint splints in the small attached hook-and-loop pockets (Fig. 3). Das „T“ (TOP) befindet sich auf dem Oberschenkel.
- Now stick the splints onto the side of the brace. Ensure that the brace joints are in the correct position: The middle of the brace joints should be level with the top of the kneecap (Fig. 4). The joints must be positioned behind the centre of the side of the leg and must run parallel to each other. The joints must be adapted to the patient's individual leg shape by limiting them.

- Wrap both straps around the leg and close the straps on the brace (Fig. 5).
- If the straps are too short, simply fasten both of the additional strap extensions on the hook-and-loop areas at the ends of the straps for the brace (Fig. 6). You can now wrap the straps around the leg as usual and fasten them on the brace.

---

### Care instructions

Please close the hook and loop fasteners before washing and remove the splint.

Switch the valve to “OFF”. We recommend using a laundry bag. Soap residues, lotions and ointments can cause skin irritation and material wear. The cooling compress can be cleaned carefully with water, a moistened towel or a sponge.

- Wash the product by hand, preferably using a medi clean detergent, or in delicate cycle at 30°C using a mild detergent without fabric conditioners.
- Do not bleach.
- Leave to dry naturally.
- Do not iron.
- Do not dry clean.

---

### Storage instructions

Keep the product in a dry place and do not expose to direct sunlight.

---

### Material composition

Knee sleeve: Cotton, nylon

Blow-up inner cushion: Polyester, nylon

Air pump: Silicone, polyethene, synthetic rubber

Cooling compress: Propylene glycol,

sodium carboxymethyl cellulose, water

Joints: Aluminium, polyester

### Liability

The manufacturer's liability will become void if the product is not used as intended. Please also refer to the corresponding safety information and instructions in this manual.

### Disposal

The product can be disposed of in the domestic waste.

In the event of any complaints regarding the product such as damage to the fabric or a fault in the fit, please report to your specialist medical retailer directly. Only serious incidents which could lead to a significant deterioration in health or to death are to be reported to the manufacturer or the relevant authorities in the EU member state. The criteria for serious incidents are defined in Section 2, No. 65 of the Regulation (EU) 2017/745 (MDR).

### Note on additional functions of the brace

The brace can also be used with a cooling compress and/or with compression with the aid of an additional air pump, dependent on the instructions of the doctor.

### Warning!

Please consult the doctor treating you before applying compression and/or cold therapy. This is mandatory. It may be that this treatment is contraindicated for you.

If the cooling compress is defective for any reason, if it leaks and the gel comes into contact with the skin or eyes, wash the affected areas thoroughly with water straight away. Never try to heat the cooling compress in any way (e.g. water baths, microwaving, ovens etc. are all prohibited).

**Contraindications for cooling and compression therapy**

If you have a vascular disease, circulation disorder, paraesthesia or tissue inflammation, these therapies may only be introduced once they have been discussed with your doctor.

**Additional function: Cooling**

- To use the cooling function, leave the included cooling compress in the freezer briefly (max. 30 minutes) or simply in the fridge. An excessively long cooling period in the fridge is not recommended as the temperature could then drop too low and injuries or frostbite could occur.
- After successful cooling, the cooling compress (Fig. 7) can simply be fastened to the black inner side of the brace (the blue inner side of the cooling compress is left facing you).

**Warning!**

Ensure that an underlay (thin towel or piece of clothing) is always present between the skin and the cooling compress. The length of time and frequency of the cooling phases must always be specified by the treating doctor and should not last longer than a max. of 30 minutes.

**Additional function: Compression**

- To apply compression, attach the tube to the valve (upper right) (Fig. 8).
- Switch the valve to "ON" (Fig. 9).
- Pump carefully until you have achieved the desired level of compression (Fig. 10). After filling has been completed, switch the valve to "OFF".
- To let some air out, you can press the button on the pump (Fig. 11).
- If you want to let all of the air out, remove the tube and set the valve setting to "ON". Now the air can escape again.

**Warning!**

Do not overfill the air chamber. Excessive pressure can cut off the blood vessels and can lead to decreased circulation of blood and other fluids.

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